

# SCALEX™ 880 PLUS

## PIEZO ULTRASONIC SCALER



- State-of-the-art piezoelectric technology is self-tuning and load-adjusting for smooth and effortless operation.
- Special TURBO function helps remove especially hard calculi by increasing power relative to the hardness.
- Streamlined compact size and slim lightweight handpiece.
- Innovative water-cooled circuitry designed for greater stability during operation.
- Extremely gentle on the patient.
- Available selection of a wide variety of high quality tips to correspond with specific clinical applications.
- Unit is available for 100-240V 50/60Hz operation.
- Three tips included.

ISO 9001  
ISO 13485



**DENTAMERICA®**

SCALEX™ 880 PLUS

# SCALEX™ 880 PLUS

## PIEZO ULTRASONIC SCALER

### ※PIEZO SCALER TIPS



Reorder No. 88006  
Perio



Reorder No. 88004  
Universal B



Reorder No. 88005  
Chisel



Reorder No. 88003  
Universal A




Reorder No. 88002  
Endo Chuck

### SPECIFICATIONS:

- INPUT: 100V-240V 50/60Hz 0.8A
- OUTPUT: DC 24V 1.5A
- FREQUENCY: 27,000 to 32,000 Hz/sec
- DIMENSIONS: W140xD140xH50mm
- WEIGHT: 1,000g

### DENTAMERICA®

P.O. Box 3200  
Industry, CA 91744 U.S.A.  
Phone: (626) 912-1388  
Fax: (626) 913-0510  
www.dentamerica.com  
www.dentamerica.com.tw

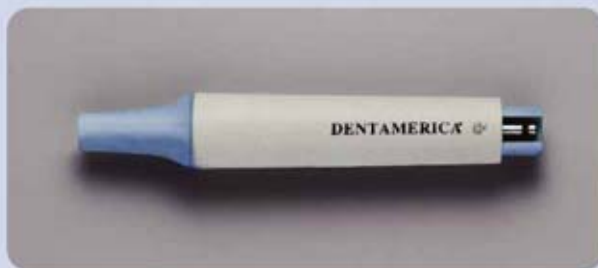

Reorder No. 88015  
Handpiece

Autoclavable handpiece with piezo ceramic crystals and titanium tip socket structure.



Reorder No. 88060  
Foot Control

Unique compact foot control designed to be activated at any angle and accessible from 360° degrees.



Reorder No. 88009  
Torque-limiting  
Tip Wrench

Special autoclavable wrench prevents damage to scaler tip when tightening.

DISTRIBUTOR

DENTAMERICA®