

LITEX® 696 TURBO

CORDLESS LED CURING LIGHT



LITEX® 696 TURBO

- Ultra high intensity over 2,000 mW/cm² cures composite in just 3 seconds/4mm.
- Equipped with both a high intensity curing mode and two stage dual-intensity ramping mode.
- Two-stage polymerization with initial lower intensity prior to higher intensity polymerization resulting in less shrinkage stress.
- Patented 360° rotatable light guide with front loaded LED provides easy access and maximum power.
- Smooth and contoured ergonomic design offers added comfort and easy cleaning.
- Large lithium battery capacity provides hundreds of curing cycles on a full charge.
- Eye protection shield included for safe operation.
- Built-in light intensity indicator.
- Low battery warning beeper.

ISO 9001
ISO 13485




GMP


DENTAMERICA®

LITEX[®] 696 TURBO

CORDLESS LED CURING LIGHT

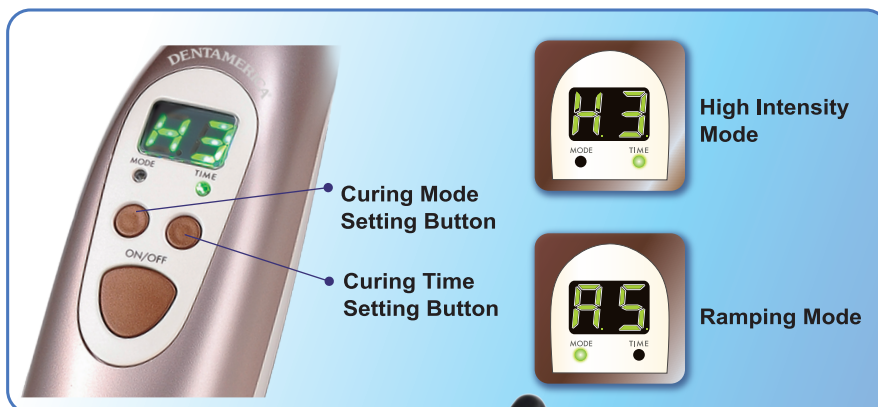

3-second cure for 4mm curing depth

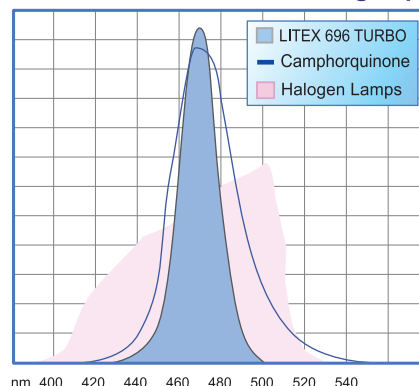





SPECIFICATIONS:

- **INPUT** : 100-240AC 0.36A
- **HANDPIECE** : 155g
L180 x W36x H40mm
- **OUTPUT** : 9V DC 1.3A
- **CHARGING UNIT** : 280g
L134x W92x H157mm



Depth of Cure

Composite tested: DENTAMERICA A3

	5 sec	3 sec	
High Intensity Mode	 5mm	 4.2mm	Depth of cure (mm)
Ramping Mode	 4.5mm	 4.1mm	




High Power Blue LED Probe



Orange Eye Protection Shield
Provides safe operation



LED Probe Cover
For LED protection

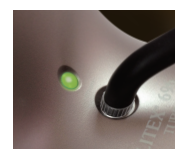

Built-in Light Intensity Indicator

DENTAMERICA[®]

18688 E. San Jose Ave.
City of Industry, CA 91748 U.S.A.
Phone: (626) 912-1388
Fax: (626) 913-0510
www.dentamerica.com
www.dentamerica.asia

DISTRIBUTOR

DENTAMERICA[®]